


BH SERIES

BIOHAZARD CLASS A HOOD WITH LAMINAR FLOW



SAFETY BEYOND ALL LIMITS

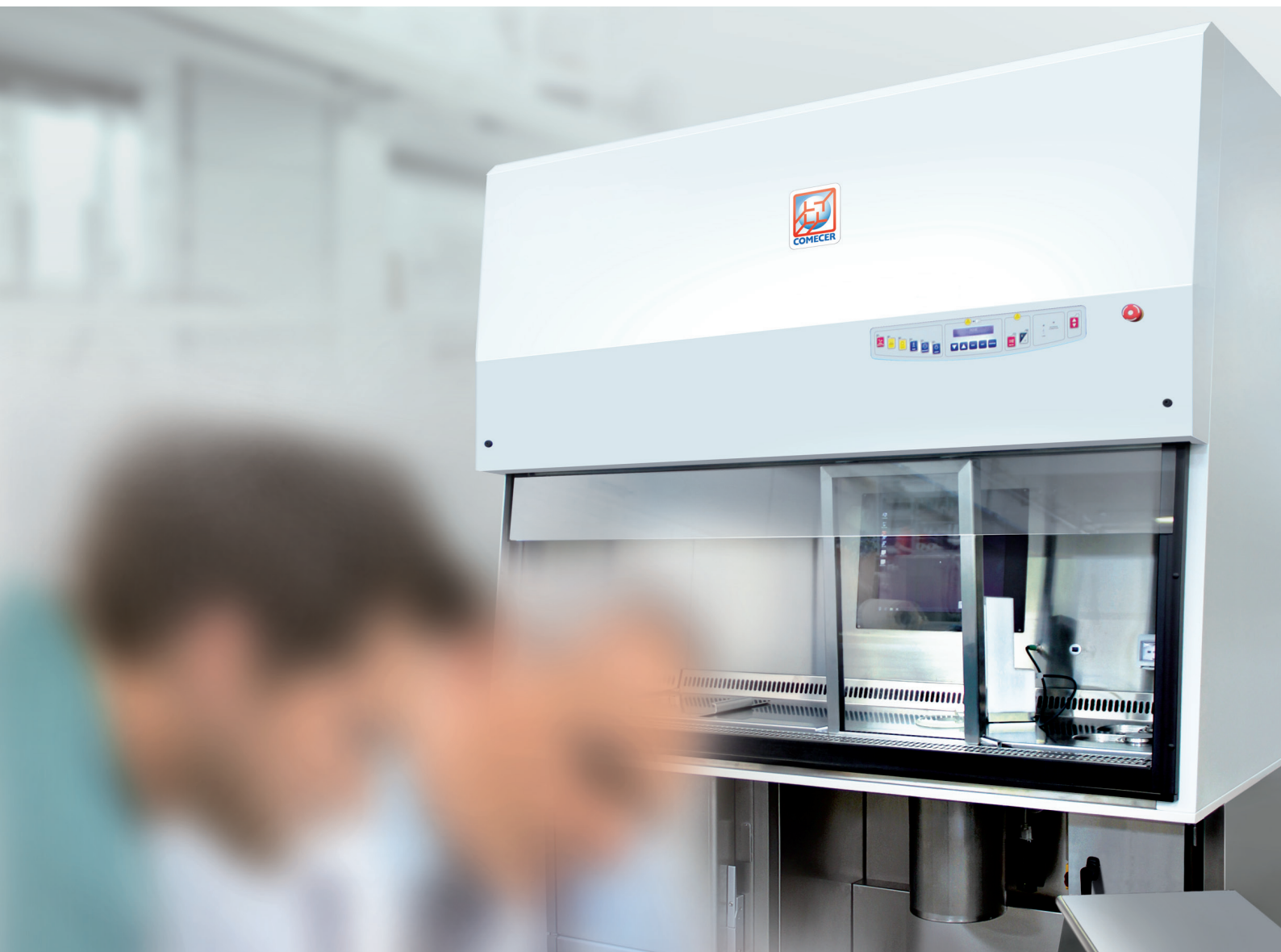

RADIOPHARMA



COMECER

an  company

BH SERIES


BIOHAZARD CLASS A HOOD WITH LAMINAR FLOW

The series of BH fume hoods is a family of workbenches with laminar flow (LAF), specifically designed to work with radiopharmaceuticals in sterile and safe conditions.

- Work area with air quality conforming to Class A "At rest" (EEC-cGMP)
- Designed specifically to work with radiopharmaceuticals
- Maximum safety for the user, the product and the environment
- Spacious work chamber for maximum ergonomics

Main technical features

- Available with Pb 5 or 10 or 30 or 50mm shielding in side panels and work surface
- Class A "At rest" (EEC-cGMP) working area with laminar flow on the entire area (downflow)
- Air inflow according to reference standards
- Filtration system to generate laminar flow in Class A area, made with HEPA H14 absolute filtering cartridge
- Air outflow filtration system made with HEPA H14 absolute filter
- Work area made of AISI 304 stainless steel
- Shielded glass window on the front side
- Dose calibrator compartment
- Shielded solid waste compartment with extraction door outside the class A area
- Shielded generator compartment to automatically lift 2 Mo99/Tc99m generators up to the work surface
- A 17" monitor in the rear side of the work area, to control the dose calibrator
- PC support and calibrator console on the front side

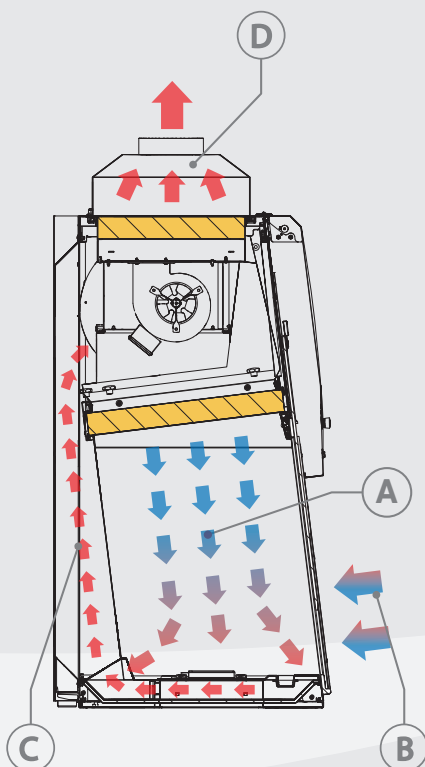

The BH Series is designed and built according to the international standards for microbiological safety cabinets as defined by EN12469 (Class II).

The pressurised air pushed by the fan passes through the absolute filter and falls, with laminar flow, into the work area (A). Here, it is mixed with the outdoor air, which enters through the front opening (B), to then be drawn through the perforated work surface, inside the recirculation channel, at the back of the work area (C).

Part of this air is expelled through the dedicated HEPA filter (D), allowing the air to enter from outside through the front barrier. This way, the operator and the surrounding environment are guaranteed to be protected from hazardous materials handled in the work area.



Working chamber

Shielded mobile protection

Generator compartment

Dose calibrator

Solid waste compartment



Work chamber - Spacious work chamber made of AISI 304 stainless steel for maximum ergonomics, equipped with a motorised safety 8 mm thick glass window. A 17" monitor integrated in the rear side of the work area, used to control the dose calibrator*.



Sliding Lead glass window* to protect the operator from radiation.



Generator compartment* - The generator compartment is separated from the work surface by a hinged door.

Each generator is handled from its own "ascent" button. The generator compartment is also equipped with a motorised shielded damper that makes sure the unused generator is perfectly shielded.



Solid waste compartment* - The waste is stored in the compartment through an opening on the work surface. The waste is then removed through the shielded front door. This process follows the GMP guidelines, according to which, the waste must not be removed directly from the sterile work area.



Dose calibrator* - The ionisation chamber compartment is at the centre of the machine work surface. The sample holder lifter automates the lifting action of the vials/syringes from the ionisation chamber to the work surface and vice-versa.

* On the customer's request

GENERAL FEATURES

Equipment lines

The machine is offered in different lines of equipment so as to meet every requirement.

Main equipment	Models	BH4		BH5		BH6	
		Standard	Reinforced	Standard	Reinforced	Standard	Reinforced
Machine Width		1350 mm	1480 mm	1655 mm	1785 mm	1960 mm	2090 mm
Protective window		S	S	S	S	S	S
Liquids collection tray and pre-filter		S	S	S	S	S	S
Work area shielding (mm of Pb)*							
5 mm of Pb		O	-	O	-	O	-
10 mm of Pb		O	-	O	-	O	-
30 mm of Pb		-	S	-	S	-	S
50 mm of Pb		-	O	-	O	-	O
Solid waste compartment							
10 mm of Pb shielding		O	O	O	O	O	O
30 mm of Pb shielding		O	O	O	O	O	O
Dose calibrator compartment							
20 mm of Pb shielding		O	O	O	O	O	O
50 mm of Pb shielding		O	O	O	O	O	O
Generator compartment							
40 mm of Pb shielding (Light)		O	O	O	O	O	O
50 mm of Pb shielding (Ga68)		O	O	O	O	O	O
50 mm of Pb shielding (Tc99)		O	O	O	O	O	O
80 mm of Pb shielding (Tc99)		O	O	O	O	O	O
Lead shielded sliding glass window							
5 mm of Pb (equivalent) shielding		O	-	O	-	O	-
10 mm of Pb (equivalent) shielding		O	-	O	-	O	-
30 mm of Pb (equivalent) shielding		-	S	-	S	-	S
50 mm of Pb (equivalent) shielding		-	O	-	O	-	O
Dose calibrator							
VDC-606 dose calibrator up to 2 Ci		O	O	O	O	O	O
VDC-606 dose calibrator up to 20 Ci		O	O	O	O	O	O
IBC Dose Calibrator up to 2 Ci		O	O	O	O	O	O
IBC Dose Calibrator up to 20 Ci		O	O	O	O	O	O

S- Standard; O- Optional;

Technical data

Structure material		AISI 304 - Scotch-Brite™	
Working chamber material		AISI 304 - Scotch-Brite™	
Hood covering material		Carbon steel treated with epoxy paints	
Lead purity	Title	Pb 98% + Sb 2%	
Weight (without shieldings)	kg	BH4	310
		BH5	365
		BH6	415
Internal dimensions	mm	BH4	1190 x 580 x 740 (w x d x h)
		BH5	1495 x 580 x 740 (w x d x h)
		BH6	1800 x 580 x 740 (w x d x h)
External dimensions (standard version)	mm	BH4	1350 x 860 x 2535 (w x d x h)
		BH5	1655 x 860 x 2535 (w x d x h)
		BH6	1960 x 860 x 2535 (w x d x h)
External dimensions (reinforced version)	mm	BH4	1480 x 927 x 2575 (w x d x h)
		BH5	1785 x 927 x 2575 (w x d x h)
		BH6	2090 x 927 x 2575 (w x d x h)



COMECER S.p.A. - Via Maestri del Lavoro, 90
48014 - Castel Bolognese (RA) - Italy
t: +39 0546 656375 - f: +39 0546 656353
comecer@comecer.com - www.comecer.com



www.comecer.com/bh-series

Vertriebspartner Deutschland



WINKGEN
medical systems

email@winkgen.de
www.winkgen.de

D-35584 Wetzlar
Eichendorffstraße 3
Fon +49 6441 381437
Fax +49 6441 381442

D-89077 Ulm
Kleiststraße 45/1
Fon +49 731 140666-0